

Angela Lost 22kg (48 lbs)*

“I visited a dietitian every month. I learned about lifestyle change, lifestyle eating changes... it’s the best thing I’ve ever done.”

The ORBERA® 2-Part Program



ORBERA®



Diet and Exercise

The ORBERA® System is a tool that is designed to give you the extra help you need to adjust to a healthier lifestyle

ORBERA® combines a clinically tested and proven medical device, with your own customised plan and support team, to effectively manage weight loss. The comprehensive two-part program starts with a soft balloon placed in your stomach for six months to encourage portion control. A team of support experts will then help guide you through a diet and exercise program. At six months, the balloon is removed and your support team will continue to guide you toward making healthy lifestyle choices. They will help you retrain your appetite, adopt new nutritional habits, and establish a reasonable exercise routine that will be essential to your long-term success.

Important ORBERA® Intra-gastric Balloon System Safety Information

The ORBERA® Intra-gastric Balloon System is a weight loss aid for adults suffering from obesity, with a body mass index (BMI) ≥ 27 kg/m² or greater, who have tried other weight loss programs, such as following supervised diets, exercise, and behavior modification programs, but who were unable to lose weight and keep it off.

To receive ORBERA® you must be willing to also follow a 12-month program, beginning with the placement of ORBERA® and continuing for 6 months after, that includes a healthy diet and exercise plan. If the diet and exercise program is not followed, you will not experience significant weight loss results, in fact, you may not experience any weight loss.

Losing weight and keeping it off is not easy, so you will be supervised throughout the program by a team of physicians, physiologists, psychologists, and nutritionists. This team will help you make and maintain major changes in your eating and exercise habits.

ORBERA® is placed for no more than 6 months. Any time that the balloon is in the stomach for longer than 6 months puts you at risk for complications, such as bowel obstruction, which can be fatal.

Some patients are ineligible to receive ORBERA®. Your doctor will ask you about your medical history and will also perform a physical examination to determine your eligibility for the device. Additionally, at the time of placement, the doctor may identify internal factors, such as stomach irritation or ulcers, which may prevent you from receiving ORBERA®.

You must not receive ORBERA® if you are pregnant, a woman planning to become pregnant within six months’ time, or breast-feeding.

Complications that may result from the use of ORBERA® include the risks associated with any endoscopic procedure and those associated with the medications and methods used in this procedure, as well as your ability to tolerate a foreign object placed in your stomach. Possible complications include: partial or complete blockage of the bowel by the balloon, insufficient or no weight loss, adverse health consequences resulting from weight loss, stomach discomfort, continuing nausea and vomiting, abdominal or back pain, acid reflux, influence on digestion of food, blockage of food entering the stomach, bacterial growth in the fluid filling the balloon which can lead to infection, injury to the lining of the digestive tract, stomach or esophagus, and balloon deflation.

Important: For full safety information please visit orbera.com/dfu or talk with your doctor.

CAUTION: Rx only.

References:


1. Directions For Use (DFU). ORBERA® Intra-gastric Balloon System (ORBERA®) Austin, TX: Apollo Endosurgery, Inc. (GRF-00346-00R01)

© 2015 Apollo Endosurgery, Inc. All rights reserved. Any third-party trademarks used herein are the property of their respective owners.

*Patient shown here is an ORBERA® patient who was treated in Australia.



MKT-00888-01 R01
Date of preparation: November 2015

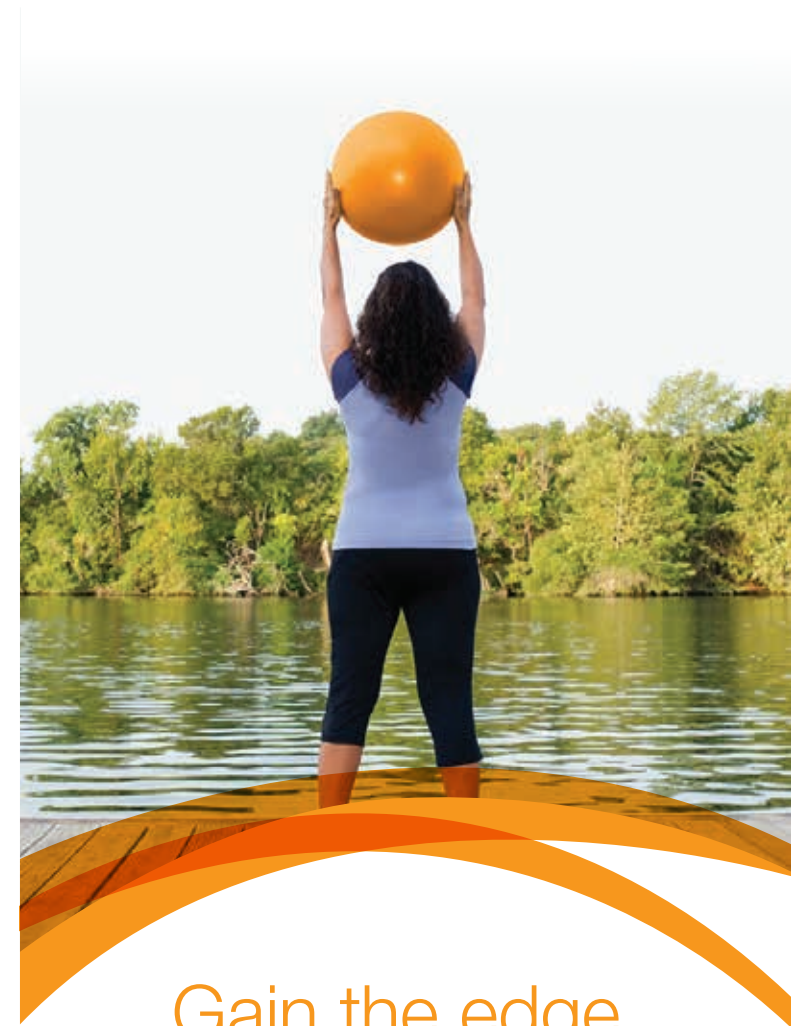

**Gain the edge.
Lose the weight.**



The ORBERA® Advantage

- After six months with the ORBERA® two-part program, people lost an average 3.1x the weight patients lost with diet and exercise alone.
- Custom support provided by a team of experts.
- Learn healthy habits to help you keep the weight off even after the balloon is removed.
- A more than 20-year history of helping thousands of people lose weight.
- #1 in global market leader. More than 220,000 devices implemented worldwide in over 80 countries.
- With the FDA approval of ORBERA®, patients now have a safe and effective treatment option.

TO LEARN MORE ABOUT ORBERA®
VISIT ORBERA.COM

What to Expect

With ORBERA®, people lost an average 3.1x the weight patients lost with diet and exercise alone. The comprehensive two-part weight loss program will help you create healthy lifestyle habits and maintain your weight loss over time.



Most weight is lost in the first three months, but because the gastric balloon stays in for six months, you will learn the right amount of food to eat to maintain your new weight.



After six months with the ORBERA® two-part program, people lost an average of 3.1x the weight than diet and exercise alone.

20-Years Global Clinical Experience

ORBERA® is a clinically tested and proven weight loss solution that has helped people in over 80 countries. No fad programs. No surgery. Just results.



More than 220,000 devices distributed worldwide.



Assisting weight loss for over 20 years.



CE marked.

Support to Meet Your Long-Term Weight Loss Goals

Once the balloon is in place, you'll receive your own customised diet and exercise plan and coaching from a support team of experts to help keep you motivated, provide advice, coordinate your program, and help you work through weight loss barriers.



ORBERA® Specialist
a licensed doctor of medicine specialising in the field of weight loss.



Professional Dietitian
an expert on human nutrition and diet focusing on areas of disease prevention and weight maintenance.



Exercise Trainer
a professional specialising in creating exercise plans designed to prevent illness and restore health and well-being.

Is ORBERA® Right For Me?

While ORBERA® does not require surgery, any medical procedure requires careful thought and consideration.

To qualify for ORBERA® you will need to:

- Be an adult (over 18 years old)
- Have a Body Mass Index (BMI) greater than 27
- Be willing to participate in a medically supervised program



If you are excited to commit to this new lifestyle, visit ORBERA.COM today.